

model **BRW1000**

BARIATRIC WHEELCHAIR SCALES

Detecto
Model BRW1000
Portable Bariatric
Wheelchair Scale




Optional

- Robust 1,000 lb / 450 kg capacity
- 40-inch extra-wide, portable platform
- Two-way, easy-access, integral ramps

Detecto's model BRW1000 portable digital wheelchair scale is designed specifically for bariatric-size wheelchair patients. The BRW1000 features a low-profile platform with two-way, easy-access integral ramps for efficient wheelchair loading. The digital scale's spacious platform measures extra-large 40-inch wide by 30-inch deep dimensions to accommodate oversized bariatric wheelchairs.

The robust 1,000-lb / 450-kilogram capacity and extra-large platform size allows the BRW1000 to be used with patients in manual or power wheelchairs, dialysis patients, or walk-ons.

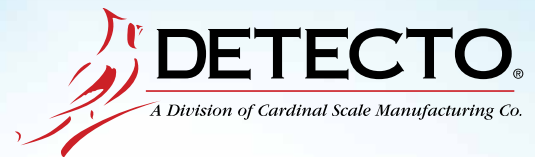

Large Bariatric Wheelchair Size Platform

This lightweight aluminum digital wheelchair scale is made for heavyweight tasks in patient healthcare and will even handle power wheelchairs.

The BRW1000's MedVue indicator is attached to the scale via a pigtail cable that extends up to 14 ft. long.



Integral Wheels combine with a lightweight frame and battery power to make the scale completely mobile.



Carrying Handles ease transport and allow the scale to be used in multiple locations.



Stainless Steel Load Cells feature self-aligning feet, so they will find the most level surface for weighing.

Weight, Height, and Body Mass Index patient measurement data may be output to EMR/EHR software through the indicator's standard RS232 serial or USB ports, or through optional wired or Wi-Fi Ethernet connectivity.

OPTIONAL WI-FI

The BRW1000 features brilliant blue LCD readouts with 7/8-in-high weight digits and up to 14-digit patient IDs that may be entered via the indicator's color-coded keypad.



The non-skid rubber mat provides traction for patients standing barefoot or wheelchair tires to grip for assurance and stability.

Optimum patient weighing "on-the-move"

- Up to 99 tares (optional) to store wheelchair weights
- Hold/Release key to retain weight after leaving the scale
- Body Mass Index calculation
- Interface with a PC or printer via multiple connectivity methods




The scale's MedVue® indicator snaps onto the platform during transport.

Clinical-Grade Construction

Integral dual wheels and two carrying handles combine with a relatively-lightweight, 75-lb frame and battery operation to allow the BRW1000 rolling mobility for versatile weighing use. Simply lift up and easily transport the BRW1000 from location to location throughout a healthcare facility wherever it's needed. The scale's four stainless steel load cells feature self-aligning feet, so when the scale comes to rest from transport mode, it will find the most level surface.

The BRW1000 is part of Detecto's Patriot line of clinical scales. From the stainless steel load cells to the MedVue® digital weight indicator, the scale is completely manufactured in Webb City, Missouri at Detecto's factory. This digital scale's quality is apparent throughout its construction.

model **BRW1000**
BARIATRIC
WHEELCHAIR
SCALES

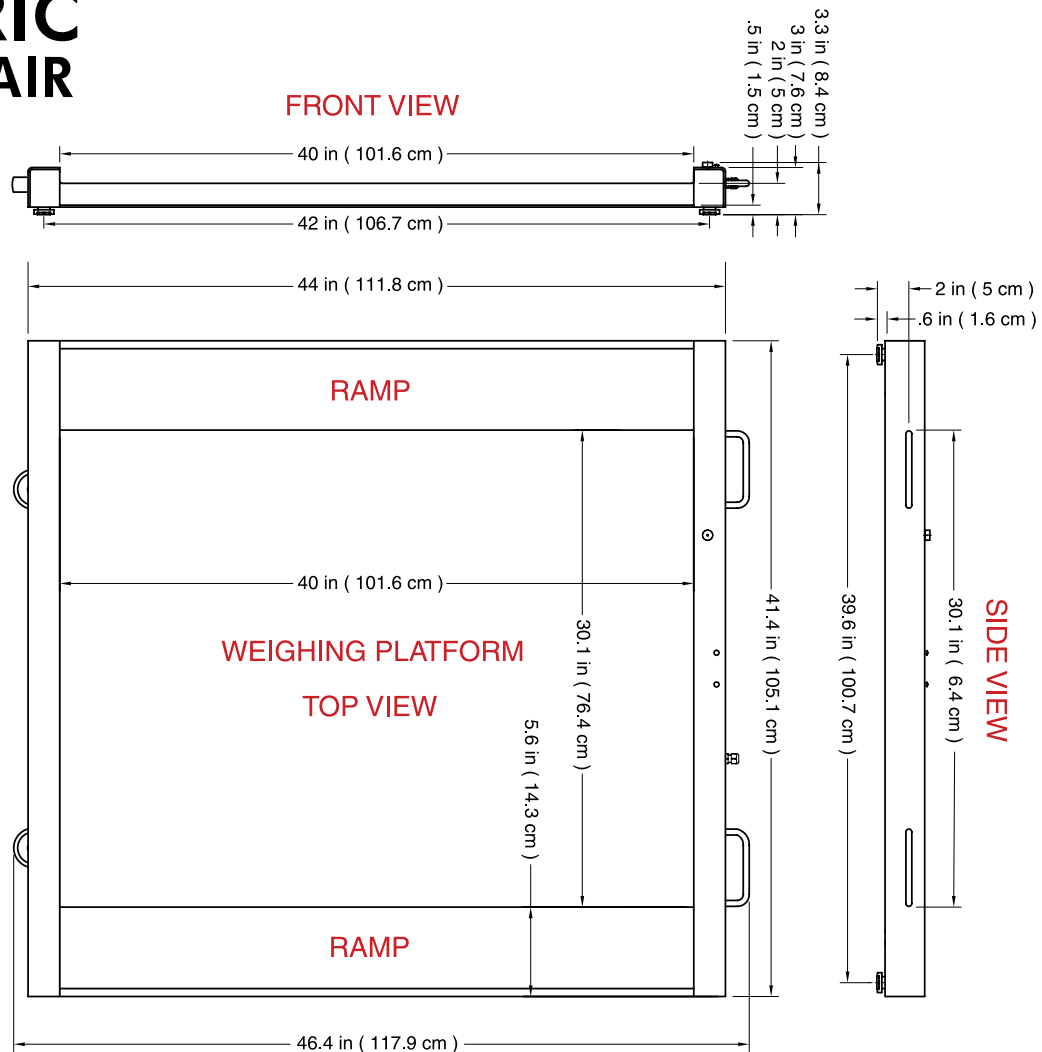

Capacity	1,000 lb x 0.2 lb / 450 kg x 0.1 kg
Display	Weight: 5-digit, 7-segment, 7/8 inch (22.2 mm) high LCD Height: 4-digit, 7-segment, 1/2 inch (12.7 mm) high LCD
Keyboard	Membrane type with 19 keys
Connectivity Ports	RS232 serial and USB (standard), wired and Wi-Fi Ethernet (optional)
Display Dimensions	6.8 inch H x 7.6 inch W x 1.8 inch D (171.4 mm W x 193.8 mm H x 44.4 mm D)
Power	6 AA cell Alkaline, Ni-Cad or NiMH batteries (not included) OR an optional 12 VDC 1.25A wall plug-in AC power adapter (Detecto part number MV1PWR)
Resolution	5,000 divisions, maximum
Operating Temperature	14° to 104° F (-10° to +40° C)
Construction	Aluminum
UPC Code	809161139303
Shipping Weight	75 lb / 34 kg

Detecto Scale reserves the right to improve, enhance, or modify features and specifications without prior notice.



DETECTO®
 A Division of Cardinal Scale Manufacturing Co.
 203 E. Daugherty, Webb City, MO 64870 USA
 Ph: 417-673-4631 or 1-800-641-2008 Fax: 417-673-5001
 www.Detecto.com

Sold By:

